DSC1
Compact Digital Servo Controller

User Guide



Table of Contents

Chapter 1 Introduction 1

1.1 Intended Use 1

1.2 Explanation of Safety Warnings 1

1.3 Description 1

1.4 Technical Data 2

1.4.1 Specifications 2

1.4.2 Mechanical Drawings 3

1.5 Simplified Declaration of Conformity 3

1.6 FCC Designation 3

Chapter 2 Safety 3

Chapter 3 Installation 4

3.1 Warranty Information 4

3.2 Included Components 4

3.3 Installation and Setup 4

3.3.1 Basics 4

3.3.2 Ground Loops and the DSC1 4

3.3.3 Powering the DSC1 5

Chapter 4 Operation 5

4.1 Touchscreen 5

4.1.1 Launching the Touchscreen Interface 5

4.1.2 Touchscreen Operation in SERVO Mode 5

4.1.3 Touchscreen Operation in RAMP Mode 7

4.1.4 Touchscreen Operation in PEAK mode 7

4.1.5 Settings Screen 9

4.2 Software 9

4.2.1 Launching the Software 9

4.2.2 Servo Software Tab 10

4.2.3 Ramp Software Tab 12

4.2.4 Peak Software Tab 13

4.2.5 Firmware Update Tab 14

4.2.6 Saved Data 15

Chapter 5 Theory of Operation 16

5.1 PID Servo Control 16

5.2 PID Theory 17

5.2.1 PID Theory for a Continuous Servo Controller 17

5.3 Bumpless Integration 18

5.4 Discrete-Time Servo Controllers 18

5.4.1 Data Format 18

5.4.2 Control Law in the Time Domain 18

5.4.3 Approximating the Integral and the Derivative 19

5.5 Servo Mode: PID Tuning 19

5.6 Ramp Mode 20
5.7 Peak Lock Mode 20
Chapter 6 Maintenance and Cleaning 23
Chapter 7 Troubleshooting and Repair 23
Chapter 8 Disposal 23
Chapter 9 Thorlabs Worldwide Contacts 24

Chapter 1 Introduction

1.1 Intended Use

The DSC1 compact digital servo controller is designed for general laboratory use in research and industry. The DSC1 controller measures a voltage, computes a feedback signal according to the user selected control algorithm, and outputs a voltage. The product may only be used in accordance with the instructions described in this manual. Any other use will invalidate the warranty. Any attempt to reprogram, disassemble binary codes, or otherwise alter the factory machine instructions in a DSC1, without Thorlabs consent, will invalidate the warranty.

Thorlabs recommends using the DSC1 with Thorlabs' photodetectors and actuators. Examples of Thorlabs actuators that are well suited to use with the DSC1 are Thorlabs' piezo amplifiers, laser diode drivers, thermoelectric cooler (TEC) controllers, and electro-optic modulators.

1.2 Explanation of Safety Warnings

NOTICE Indicates information considered important, but not hazard-related, such as possible damage to the product.

The CE/UKCA markings on the product are the manufacturer's declaration that the product complies with the essential requirements of the relevant European health, safety, and environmental protection legislation.



The wheelie bin symbol on the product, the accessories or packaging indicates that this device must not be treated as unsorted municipal waste but must be collected separately.

1.3 Description

Thorlabs' DSC1 Digital Servo Controller is an instrument for feedback control of electro-optical systems. The device measures an input voltage, determines an appropriate feedback voltage through one of several control algorithms, and applies this feedback to an output voltage channel. Users can choose to configure the operation of the device through either the integrated touchscreen display, a remote desktop PC graphical user interface (GUI), or a remote PC software development kit (SDK).

The servo controller samples voltage data with 16-bit resolution through a coaxial SMB input port at 1 MHz. To provide more accurate voltage measurements, arithmetic circuitry within the device averages every two samples for an effective sample rate of 500 kHz. The digitized data is processed by a microprocessor at high-speed using digital signal processing (DSP) techniques. The user may choose between SERVO and PEAK control algorithms. Alternatively, the user may test a systems response to DC voltage to determine the servo setpoint with the RAMP operating mode, which outputs a sawtooth wave synchronous with the input. The input channel has a typical bandwidth of 120 kHz. The output channel has a typical bandwidth of 100 kHz. The -180 degree phase lag of the input-to-output voltage transfer function of this servo controller is typically 60 kHz.

1.4 Technical Data

1.4.1 Specifications

Operating Specifications	
System Bandwidth	DC to 100 kHz
Input to Output -180 Degree Frequency[a]	>58 kHz (60 kHz Typical)
Nominal Input Sampling Resolution	16 Bit
Nominal Output Resolution	12 Bit
Maximum Input Voltage	±4 V
Maximum Output Voltage[b]	±4 V
Maximum Input Current	100 mA
Average Noise Floor	-120 dB V^2/Hz
Peak Noise Floor	-105 dB V^2/Hz
Input RMS Noise[c]	0.3 mV
Input Sampling Frequency	1 MHz
PID Update Frequency[d]	500 kHz
Peak Lock Modulation Frequency Range	100 Hz - 100 kHz in 100 Hz Steps
Input Termination	1 MΩ
Output Impedance[b]	220 Ω

a. This is the frequency at which the output reaches a -180 degree phase shift relative to the input.
b. The output is designed for connection to high-Z (>100 kΩ) devices. Connecting devices with lower input termination, R_{dev}, will reduce the output voltage range by $R_{dev}/(R_{dev} + 220\ \Omega)$ (e.g., a device with 1 kΩ termination will give 82% of the nominal output voltage range).
c. The integration bandwidth is 100 Hz - 250 kHz.
d. A low-pass filter reduces digitization artifacts in output control variable, resulting in an output bandwidth of 100 kHz.

Electrical Requirements	
Supply Voltage	4.75 - 5.25 V DC
Supply Current	750 mA (Max)
Temperature Range[a]	0 °C to 70 °C

a. Temperature range over which the device may be operated without damage. Optimal operation occurs when near room temperature.

System Requirements	
Operating System	Windows 10® (Recommended) or 11, 64 Bit Required
Memory (RAM)	4 GB Minimum, 8 GB Recommended
Storage	300 MB (Min) of Available Disk Space
Interface	USB 2.0
Minimum Screen Resolution	1200 x 800 Pixels

1.4.2 Mechanical Drawings

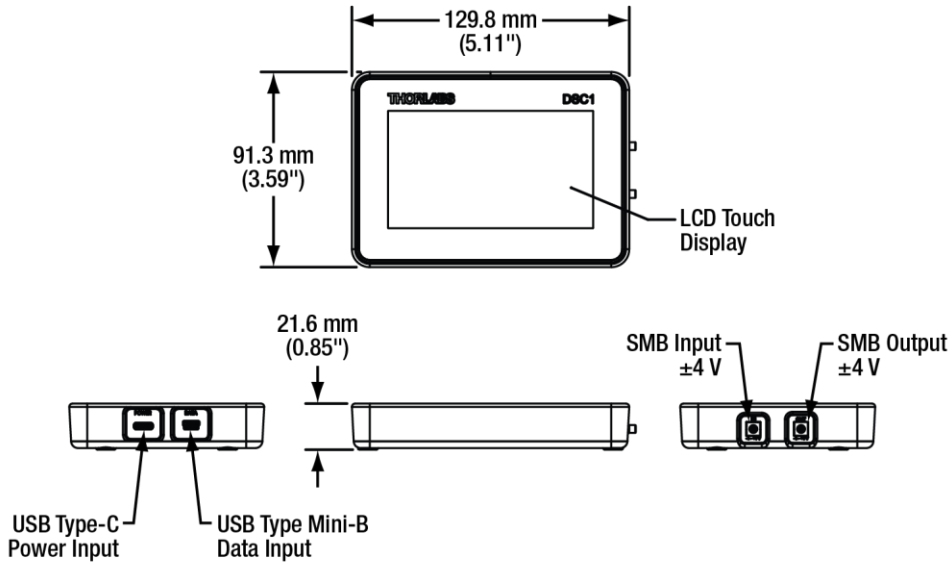


Figure 1 Mechanical Drawing

1.5 Simplified Declaration of Conformity

The full text of the EU declaration of conformity is available at the following internet address: **https://Thorlabs.com/newgrouppage9.cfm?objectgroup_id=16794**

1.6 FCC Designation

Note: This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.

Chapter 2 Safety

NOTICE

This instrument should be kept clear of environments where liquid spills or condensing moisture are likely. It is not water resistant. To avoid damage to the instrument, do not expose it to spray, liquids, or solvents.

Chapter 3 Installation

3.1 Warranty Information

This precision device is only serviceable if returned and properly packed into the complete original packaging including the complete shipment plus the cardboard insert that holds the enclosed devices. If necessary, ask for replacement packaging. Refer servicing to qualified personnel.

3.2 Included Components

The DSC1 Compact Digital Servo Controller is delivered with the following components:

- DSC1 Digital Servo Controller
- Quick Start Card
- USB-AB-72 USB 2.0 Type-A to Mini-B Data Cable, 72" (1.83 m) Long
- USB Type-A to USB Type-C Power Cable, 1 m (39") Long
- PAA248 SMB to BNC Coaxial Cable, 48" (1.22 m) Long (Qty. 2)

3.3 Installation and Setup

3.3.1 Basics

Users may configure the device with a computer using the USB interface or through the integrated touchscreen. In either case, power must be provided through the 5V USB-C connection. When using the desktop GUI, the servo controller must be connected with a USB 2.0 cable (included) from data port of the device to a PC with the Digital Servo Controller software installed.

3.3.2 Ground Loops and the DSC1

The DSC1 controller includes internal circuitry to limit the likelihood of ground loops occurring. Thorlabs suggests using either the transformer isolated DS5 regulated power supply or the CPS1 external battery pack. With either the DS5 or CPS1 power supplies, the signal ground within the DSC1 floats with respect to the earth ground of a wall outlet. The only connections to the device that are common to this signal ground are the signal ground pin of the USB-C power connector and the outer, return path on the output SMB coaxial cable. The USB data connection is isolated. The input signal has a ground-loop break resistor between the signal return path and the signal ground within the instrument which typically prevents ground loop interference. Importantly, there are no two direct paths to the device signal ground, minimizing the occurrence of ground loops. To further mitigate the risk of ground-loop interference, Thorlabs suggests the following best practices:

- Keep all power and signal cables to the device short.
- Use either a battery (CPS1) or transformer isolated (DS5) power supply with the DSC1. This ensures a floating device signal ground.
- Do not connect other instruments' signal return paths to one another.
    - A common example is a typical benchtop oscilloscope; most often the outer shells of the BNC input connections are directly connected to earth ground. Multiple ground clips connected to the same ground node in an experiment can cause a ground loop.

Although the DSC1 controller is unlikely to cause a ground loop in of itself, other instruments in a user's lab may not have ground loop isolation and thus could be a source of ground loops.

3.3.3 Powering the DSC1

The DSC1 Digital Servo Controller requires 5 V power through the USB-C at up to 0.75 A peak current and 0.55 A in typical operation. Thorlabs offers two compatible power supplies: Item #s CPS1 and DS5. In applications where noise sensitivity is less constrained or where runtimes of greater than 8 hours are required, the DS5 regulated power supply is recommended. The CPS1 battery power supply is recommended when optimal noise performance is desired. With the CPS1 fully charged and in good health, the DSC1 can operate for 8 hours or more without recharging.

Chapter 4 Operation

4.1 Touchscreen

4.1.1 Launching the Touchscreen Interface

After being connected to power and a brief, less than one second warmup, the DSC1 will illuminate the integrated touchscreen display, and the screen will respond to inputs.

In servo, peak, and ramp modes, the screen will contain a plot displaying the live data selected. The data available to be selected for display changes by mode and is described in subsequent sections.

- The color of the line in the chart matches the selected variables' highlight color - for example, in all three modes the golden color of the PV line in the plot matches the golden color highlighting the PV numeric value.
- By default, the plot is in full scale mode, displaying the entire Y-range of -4 V to +4 V. By clicking on the scale icon ( ) the chart will toggle between full scale and auto-scale mode. In auto-scale mode the plot will contract the Y-axis to the range of the current data.

4.1.2 Touchscreen Operation in SERVO Mode

The SERVO mode implements a PID controller.

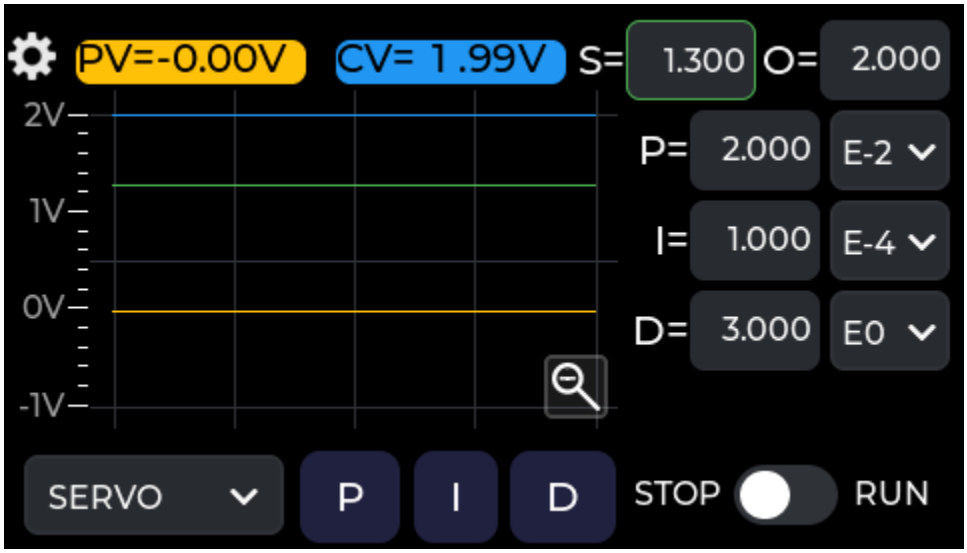


Figure 2 Touchscreen Display in the Servo Operating Mode with the PID Controller Enabled in PI Control Mode

- The **PV** (process variable) numeric value shows the average voltage of the input signal in a gold color. Pressing on the PV indicator displays or hides the process variable display in the plot.
- The **CV** (control variable) numeric value shows the current output voltage from the DSC1 in a blue color. Pressing on the CV indicator displays or hides the control variable display in the plot.
- The **S** (setpoint) control sets the setpoint of the servo loop in volts. +4 V is the maximum and –4 V is the minimum allowable, but if the setpoint is set to a value outside the current CV limits (see section 4.1.5 below) the control variable will be clipped at the maximum or minimum value and the setpoint field will be highlighted with a red background
- The **O** (offset) control sets the DC offset of the servo loop in volts. +4 V is the maximum and -4 V is the minimum allowable, but if the offset is set to a value outside the current CV limits (see section 4.1.5 below) the control variable will be clipped at the maximum or minimum value and the offset field will be highlighted with a red background
- The **P** (proportional) control sets the proportional gain coefficient. This may be a positive or negative value between 10^{-5} and 10^{5}, notated in engineering notation. See Chapter 5 for details on the absolute scaling of this coefficient.
- The **I** (integral) control sets the integral gain coefficient. This may be a positive or negative value between 10^{-5} and 10^{5}, notated in engineering notation. See Chapter 5 for details on the absolute scaling of this coefficient.
- The **D** (derivative) control sets the derivative gain coefficient. This may be a positive or negative value between 10^{-5} and 10^{5}, notated in engineering notation. See Chapter 5 for details on the absolute scaling of this coefficient.
- The **STOP-RUN** toggle disables and enables the servo loop.
- The **P**, **I**, and **D** buttons enable (illuminated green) and disable (dark blue) each gain stage in the PID servo loop.
- The **SERVO** dropdown menu allows the user to select the operating mode.
- The green trace shows the current setpoint. Each point is 2 µs apart on the X-axis.
- The golden trace shows the current measured PV. Each point is 2 µs apart on the X-axis.
- The Settings button  takes you to a Settings and Information page described in Section 4.1.5.

4.1.3 Touchscreen Operation in RAMP Mode

The RAMP mode outputs a sawtooth wave with user configurable amplitude and offset.

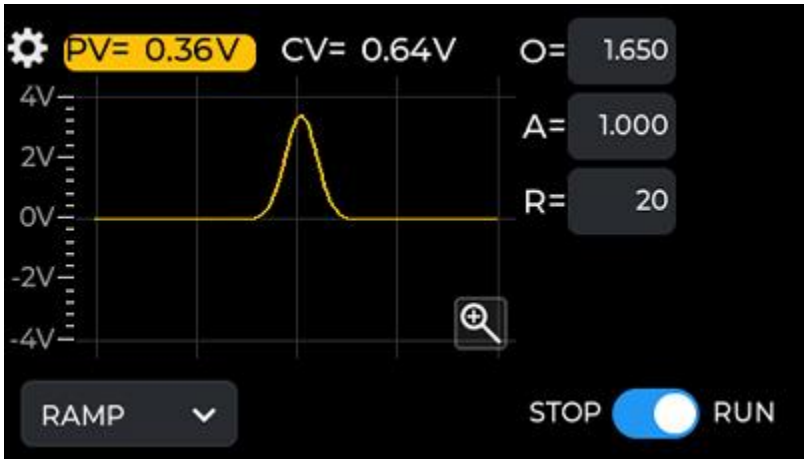


Figure 3 Touchscreen Display in the Ramp Operating Mode

- The **PV** (process variable) numeric value shows the average voltage of the input signal in volts. Pressing on the PV indicator displays or hides the process variable display in the plot.
- The **CV** (control variable) numeric value shows the current output voltage applied by the device. Pressing on the CV indicator displays or hides the control variable display in the plot.
- The **O** (offset) control sets the DC offset of the ramp output in volts. +4 V is the maximum and -4 V is the minimum allowable, but if the offset is set to a value outside the current CV limits (see section 4.1.5 below) the control variable will be clipped at the maximum or minimum value and the offset field will be highlighted with a red background.
- The **A** (amplitude) control sets the amplitude of the ramp output in volts. The sawtooth signal generated will cycle from (offset – amplitude) to (offset + amplitude). +4 V is the maximum and -4 V is the minimum allowable, but if the amplitude is set to a value outside the current CV limits (see section 4.1.5 below) the control variable will be clipped at the maximum or minimum value and the amplitude field will be highlighted with a red background.
- The **R** (ramp time) control sets the time for each ramp in milliseconds. 400 ms is the maximum and 10 ms is the minimum allowable. The **STOP-RUN** toggle disables and enables the ramp respectively. Regardless of the ramp time setting, the plot trace always shows one entire ramp cycle.
- The **RAMP** dropdown menu allows the user to select the operating mode.
- The gold-colored plot trace shows the control variable on the x-axis and the process variable on the y-axis.

4.1.4 Touchscreen Operation in PEAK mode

The PEAK mode implements an extremum seeking controller with user configurable modulation frequency, amplitude, and integration constant. Note that the modulation and demodulation are always active when the device is in PEAK mode; the run-stop toggle actives and deactivates the integral gain in the dither control loop.


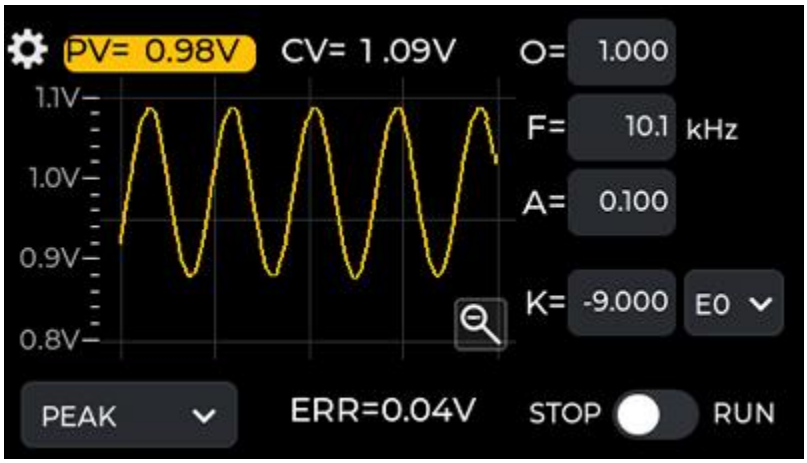


Figure 4 Touchscreen Display in the Peak Operating Mode

- The **PV** (process variable) numeric value shows the average voltage of the input signal in a gold color. Pressing on the PV indicator displays or hides the process variable display in the plot.
- The **CV** (control variable) numeric value shows the current output voltage applied by the device in a blue color. Pressing on the CV indicator displays or hides the control variable display in the plot.
- The **ERR** numeric value shows the RMS average voltage of the peak error signal in a green color. Pressing on the ERR indicator displays or hides the peak error display in the plot.
- The **F** (modulation frequency) numeric value shows the modulation frequency. The maximum modulation frequency is 100 kHz. In general, higher modulation frequencies are advisable, provided that the control actuator is responsive at that frequency.
- The **A** (amplitude) control sets the amplitude of the modulation in volts. +4 V is the maximum and -4 V is the minimum allowable, but if the offset is set to a value outside the current CV limits (see section 4.1.5 below) the control variable will be clipped at the maximum or minimum value and the amplitude field will be highlighted with a red background.
- The **K** (peak lock integral coefficient) control sets the integration constant of the controller, with units of V / s, notated in engineering notation. If the user is unsure of how to configure this value, typically starting with a value around 1 is advisable.
- The **STOP-RUN** toggle disables and enables the servo loop respectively.
- The **PEAK** dropdown menu allows the user to select the operating mode.
- The golden trace shows the process variable synchronized with the control variable. Each point is spaced 2 µs apart on the X-axis.

4.1.5 Settings Screen

Pressing on the “gear” icon in the top left brings up the settings screen:

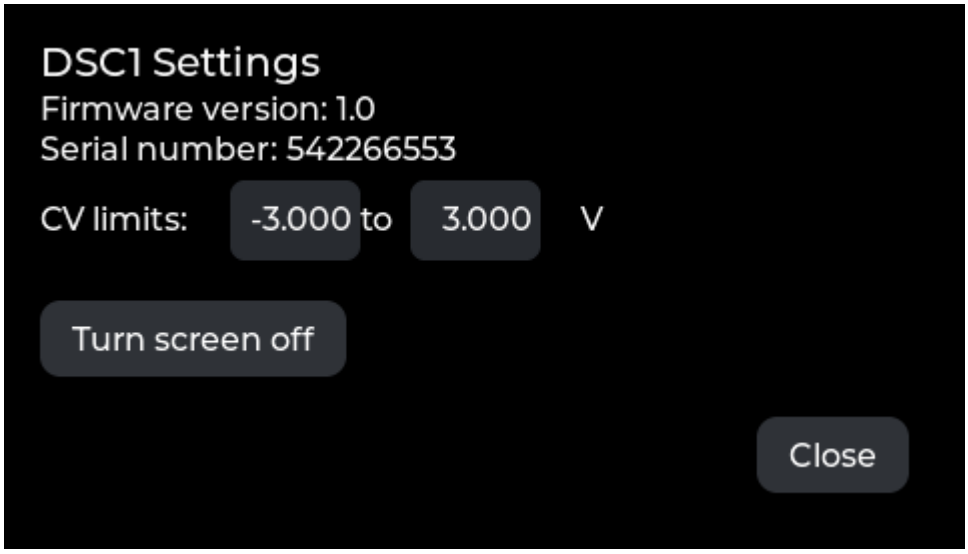

On this screen you can view the firmware version and serial number, set limits for the control variable, and turn the screen off for use in light-sensitive environments. Clicking anywhere on the screen turns the screen back on.

4.2 Software

The digital servo controller software is designed to both allow for control over basic functionality via a computer interface and provides an expanded set of analysis tools for using the controller. For example, the GUI includes a plot that can display the input voltage in frequency domain. Additionally, data can be exported as a .csv file.

This software allows for use of the device in the servo, peak, or ramp modes with control over all parameters and settings. The system response may be viewed as the input voltage, error signal, or both, either in the time domain or frequency domain representations. Please see the manual for more information.

4.2.1 Launching the Software

After launching the software, click “Connect” to list available DSC devices. Multiple DSC devices may be controlled at a time.

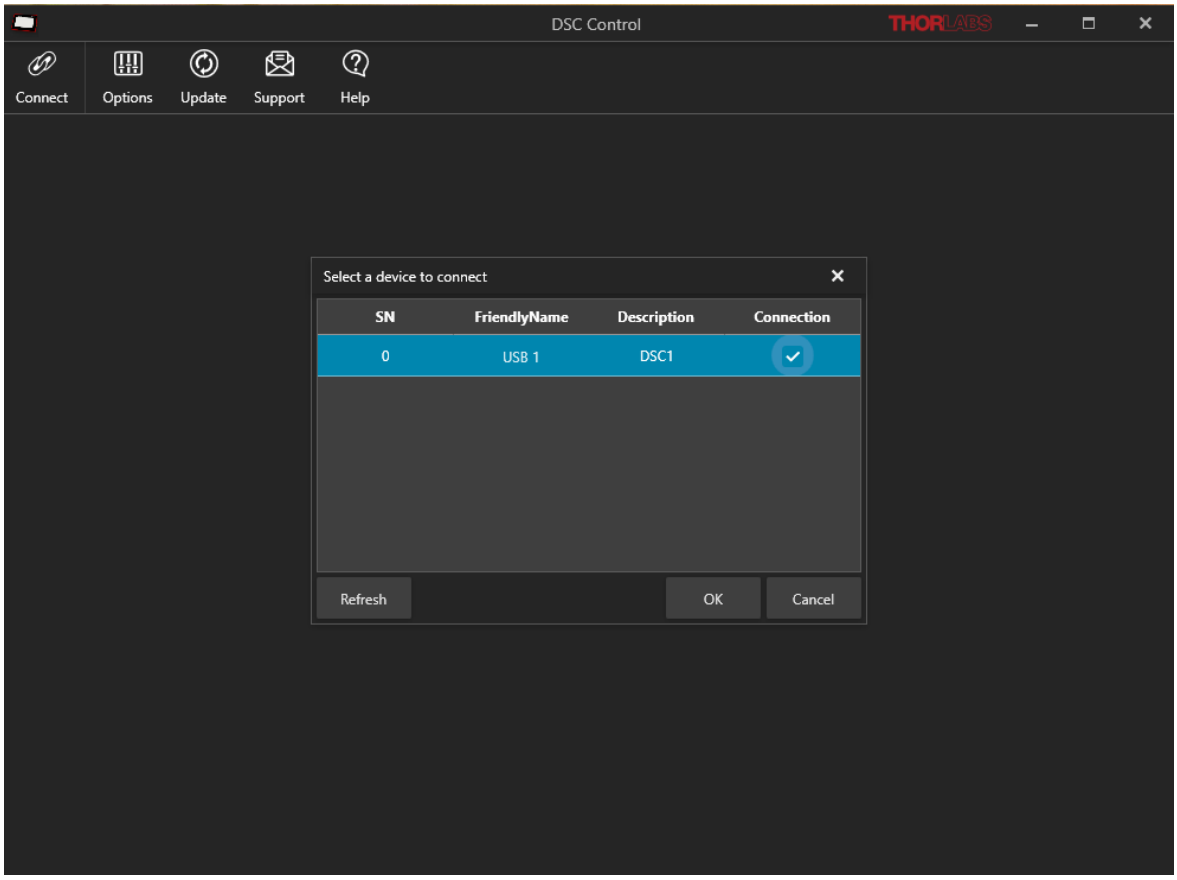


Figure 5 Launch Screen for the DSCX Client Software

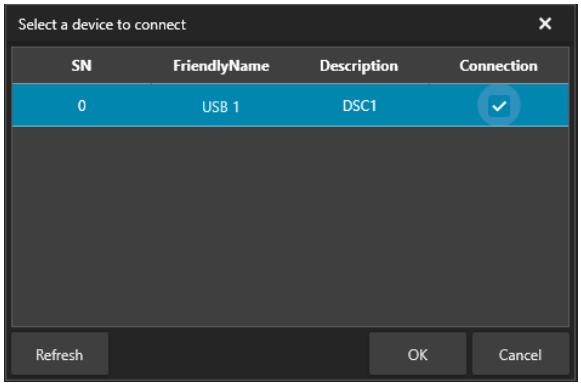


Figure 6 Device selection window - click OK to connect to the selected device.

4.2.2 Servo Software Tab

The Servo tab allows a user to operate the device in servo mode with additional controls and displays beyond those provided by the embedded touchscreen user interface on the device itself. On this tab, either time or frequency domain representations of the process variable are available. The system response may be viewed as either the process variable, error signal, or both. The error signal is the difference between the process variable and the setpoint. Using control analysis techniques, the impulse response, frequency response, and phase response of the device can be predicted, provided certain assumptions about the system's behavior and the gain coefficients are made. This data is displayed on the servo control tab so that users can preemptively configure their system, before beginning control experiments.

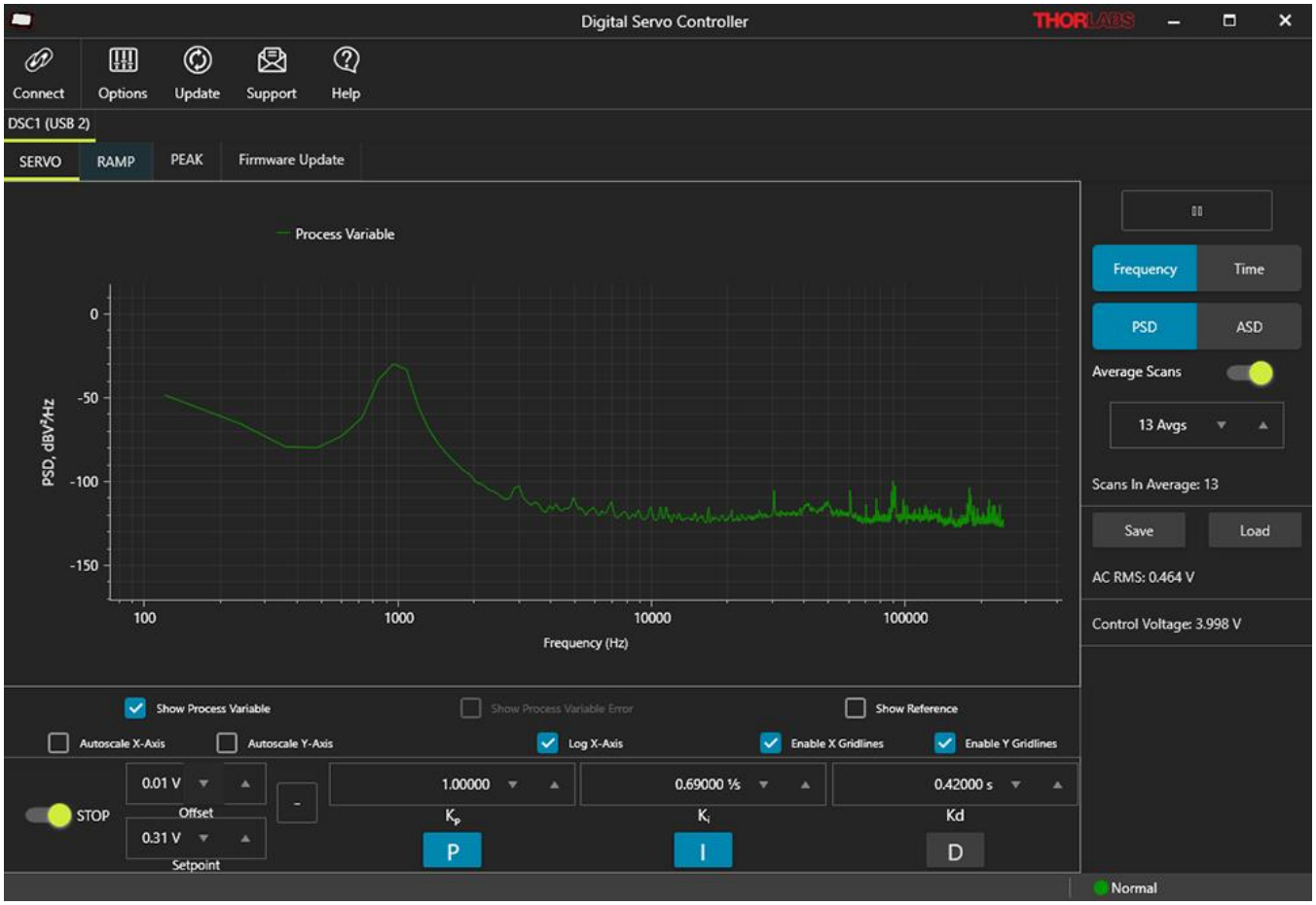


Figure 7 Software Interface in Servo Mode with the Frequency-Domain Display

- **Enable X Gridlines**: Checking the box enables X gridlines.
- **Enable Y Gridlines:** Checking the box enables Y gridlines.
- **Run / Pause Button:** Pressing this button commences / stops the update of graphical information on the display.
- **Frequency / Time Toggle:** Switches between frequency-domain and time-domain plotting.
- **PSD / ASD Toggle:** Switches between power spectral density and amplitude spectral density vertical axes.
- **Average Scans:** Toggling this switch enables and disables averaging in the frequency domain.
- **Scans In Average**: This numeric control determines the number of scans to be averaged. The minimum is 1 scan and the maximum is 100 scans. The up and down arrows on a keyboard increase and decrease the number of scans in the average. Similarly, the up and down buttons adjacent to the control increase and decrease the number of scans in the average.
- **Load**: Pressing this button in the Reference Spectrum panel allows a user to select a reference spectrum saved on the client PC.
- **Save:** Pressing this button in the Reference Spectrum panel allows a user to save the currently displayed frequency data to their PC. After clicking this button, a save file dialog will allow the user to choose the

storage location and enter the file name for their data. The data saves as a Commas Separated Value (CSV).

- **Show Reference:** Checking this box enables display of the last selected reference spectrum.
- **Autoscale Y-Axis**: Checking the box enables automatic setting of the Y Axis display limits.
- **Autoscale X-Axis:** Checking the box enables automatic setting of the X Axis display limits.
- **Log X-Axis:** Checking the box toggles between a logarithmic and linear X axis display.
- **Run PID:** Enabling this toggle enables the servo loop on the device.
- **O Numeric:** This value sets the offset voltage in volts.
- **SP Numeric:** This value sets the setpoint voltage in volts.
- **K_p Numeric:** This value sets the proportional gain.
- **K_i Numeric:** This value sets the integral gain in 1/s.
- **K_d Numeric:** This value sets the derivative gain in s.
- **P, I, D buttons:** These buttons enable the proportional, integral, and derivative gain respectively when illuminated.
- **Run / Stop Toggle:** Toggling this switch enables and disables the control.

The user may also use the mouse to change the extent of the displayed information:

- The mouse wheel zooms the plot in and out towards the current position of the mouse pointer.
- SHIFT + Click changes the mouse pointer to a plus sign. Thereafter the left-mouse button will zoom in on the position of the mouse pointer by a factor of 3. The user may also drag and select a region of the plot to zoom to fit.
- ALT + Click changes the mouse pointer to a minus sign. Thereafter the left-mouse button will zoom out from the position of the mouse pointer by a factor of 3.
- Spread and pinch gestures on a mouse pad or touch screen will zoom in and out of the plot respectively.
- After scrolling, clicking the left-mouse button will allow the user to pan by dragging the mouse.

4.2.3 Ramp Software Tab

The Ramp tab provides comparable functionality to the ramp tab on the embedded touchscreen display. Switching to this tab puts the connected device in ramp mode.

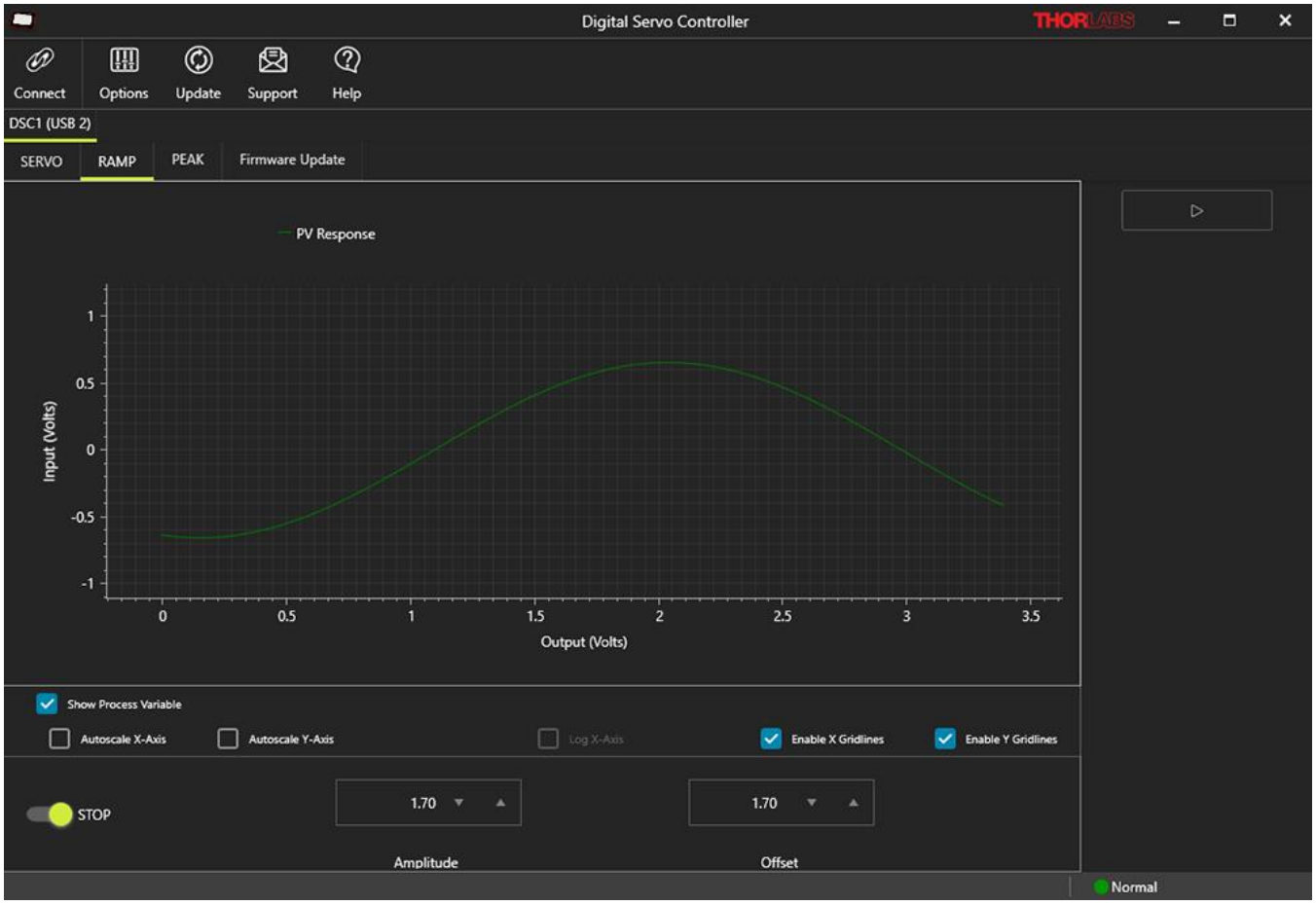


Figure 8 Software Interface in Ramp Mode.

In addition to the controls available in Servo mode, the Ramp mode adds:

- **Amplitude Numeric:** This value sets the scan amplitude in volts.
- **Offset Numeric:** This value sets the scan offset in volts.
- **Run / Stop Ramp Toggle:** Toggling this switch enables and disables the ramp.

4.2.4 Peak Software Tab

The Peak Control tab provides the same functionality as the PEAK mode on the embedded user interface, with additional visibility into the nature of the return signal from the system. Switching to this tab switches the connected device to the PEAK mode of operation.



Figure 9 Software Interface in Peak Mode with the Time-Domain Display

In addition to the controls available in Servo mode, the Peak mode adds:

- **Amplitude numeric:** This value sets the modulation amplitude in volts.
- **K numeric:** This is the peak lock integral coefficient; the value sets the integral gain constant in V/s.
- **Offset numeric:** This value sets the offset in volts.
- **Frequency numeric:** This sets the modulation frequency multiplier in increments of 100 Hz. The minimum allowable value is 100 Hz is the maximum is 100 kHz.
- **Run / Stop Peak toggle:** Toggling this switch enables and disables the integral gain. Note, whenever the device is in PEAK mode, the output modulation and error signal demodulation is active.

4.2.5 Firmware Update Tab

The Firmware Update Tab allows a user to update the firmware on the DSC1.

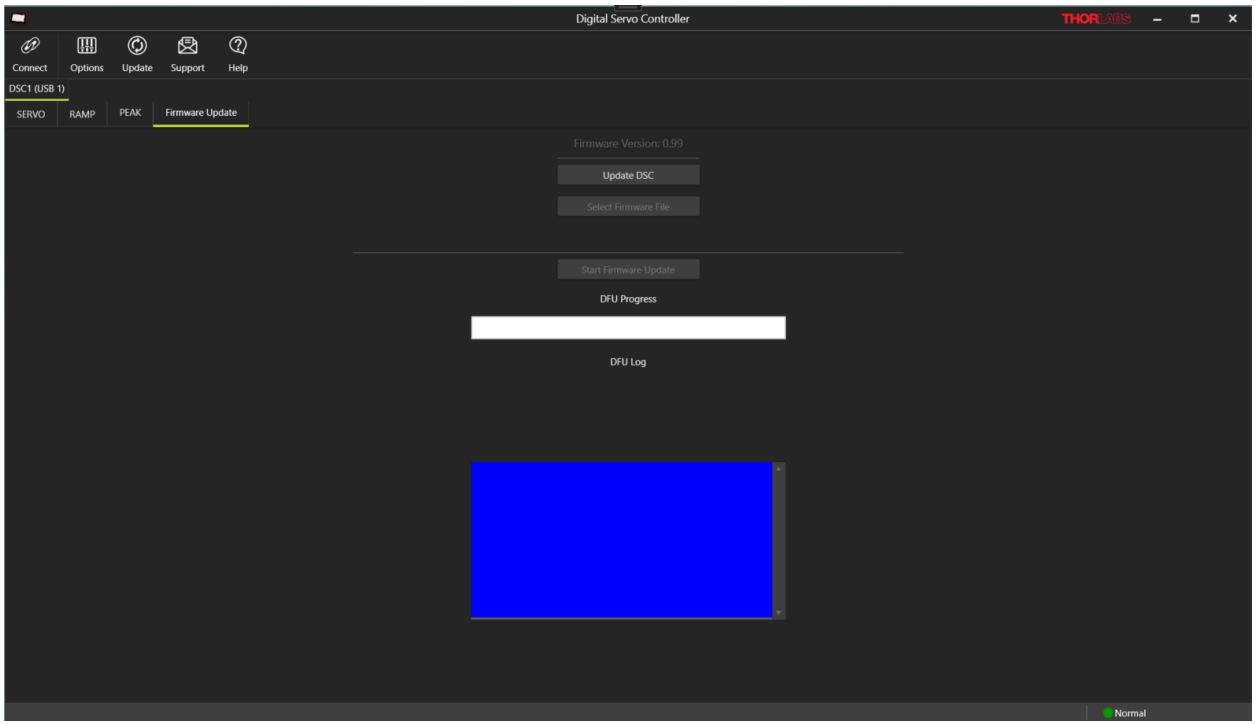


Figure 10 Software Interface in the Firmware Update Mode

To start the firmware update process, connect the DSC1 controller and check the Firmware Version shown against the latest version reported on the Thorlabs website. If a new version is available, download the TFF file.

Click on the Update DSC button to begin the process, then click on Select Firmware File and select the TFF file. Finally, click on Start Firmware Update to download the new firmware to the DSC1 controller. Do not turn off the device or close the software while the firmware update is taking place. The device should start normally with the new firmware after the update is complete. If you have any issues with the firmware update, first try power cycling the device. If problems persist, please contact **Tech Support** for assistance.

4.2.6 Saved Data

Data is saved in Comma Separated Value (CSV) format. A brief header retains pertinent data from the data being saved. If the format of this CSV is altered, the software may be unable to recover a reference spectrum. Therefore, the user is encouraged to save their data in a separate spreadsheet file if they intend to do any independent analysis.

	A	B	C	D	E
1	Timestamp	09.16.2024:09:19:11			
2	Scans an avera	1			
3	Termination	10000			
4	AC_RMS	0.000246982			
5	V_DC	-8.23E-08			
6	Frequency (Hz)	PSD (dBV^2 / Hz)	ASD (dBV / sqrt(Hz))	Sample Point(Time Domain)	Volts(Time Domain)
7	0	-130.3402866	-65.17014329	0	-0.000366211
8	122.0703	-120.6731671	-60.33658353	2	-0.000488281
9	244.1406	-119.9602799	-59.98013993	4	-0.000732422
10	366.2109	-124.711872	-62.35593601	6	-0.000610352
11	488.2813	-130.7534212	-65.37671061	8	-0.000488281
12	610.3516	-127.0521684	-63.52608422	10	-0.000244141
13	732.4219	-124.1439653	-62.07198267	12	-0.00012207
14	854.4922	-120.735596	-60.36779799	14	-0.000610352


Figure 11 Data in .csv Format Exported from the DSC1

Chapter 5 Theory of Operation

5.1 PID Servo Control

The PID circuit is often utilized as a control loop feedback controller and is very common in servo circuits. The purpose of a servo circuit is to hold the system at a predetermined value (set point) for prolonged periods of time. The PID circuit actively holds the system at the set point by generating an error signal that is the difference between the set point and the current value and modulating an output voltage to maintain the set point.

The letters making up the acronym PID correspond to Proportional (P), Integral (I), and Derivative (D), which represent the three control settings of a PID circuit. The proportional term is dependent upon the present error, the integral term is dependent upon the accumulation of past error, and the derivative term is the prediction of future error. Each of these terms are fed into a weighted sum that adjusts the output voltage of the circuit, u(t). This output is fed into the control device, its measurement is fed back into the PID loop, and the process is allowed to actively stabilize the circuit's output to reach and hold the set point value. The block diagram below illustrates the action of a PID circuit. One or more of the controls can be utilized in any servo circuit depending on what is needed to stabilize the system (i.e., P, I, PI, PD, or PID).

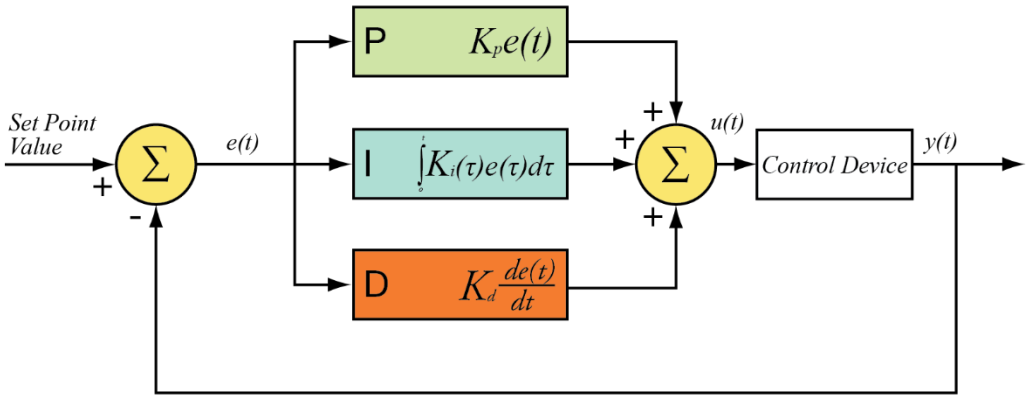

Please note that a PID circuit will not guarantee optimal control. Improper setting of the PID controls can cause the circuit to oscillate significantly and lead to instability in control. It is up to the user to properly adjust the PID parameters to ensure proper performance.

5.2 PID Theory

5.2.1 PID Theory for a Continuous Servo Controller

The output of the PID control circuit, u(t), is given as

$$u(t) = K_P\, e(t) + \int_0^t K_i(\tau)e(\tau)d\tau + K_d \frac{d}{dt}e(t)$$

Where:

K_P is the proportional gain, dimensionless

$K_i(\tau)$ is the integral gain in 1/seconds

K_d is the derivative gain in seconds

$e(t)$ is the error signal in volts

$u(t)$ is the control output in volts

From here we can define the control units mathematically and discuss each in a little more detail. Proportional control is proportional to the error signal; as such, it is a direct response to the error signal generated by the circuit:

$$P = K_P e(t).$$

Larger proportional gain results in larger changes in response to the error and thus affects the speed at which the controller can respond to changes in the system. While a high proportional gain can cause a circuit to respond swiftly, too high a value can cause oscillations about the set point value. Too low a value and the circuit cannot efficiently respond to changes in the system.

Integral control goes a step further than proportional gain, as it is proportional to not just the magnitude of the error signal but also the duration of any accumulated error.

$$I = \int_0^t K_i(\tau)e(\tau)d\tau$$

Integral control is highly effective at increasing the response time of a circuit along with eliminating the steady-state error associated with purely proportional control. In essence, integral control sums over any previously uncorrected error multiplied by K_i to produce the integral response. Thus, for even a small, sustained error, a large aggregated integral response can be realized. However, due to the fast response of integral control, high gain values can cause significant overshoot of the setpoint value and lead to oscillation and instability. Too low and the circuit will be significantly slower in responding to changes in the system.

Derivative control attempts to reduce the overshoot and ringing potential from proportional and integral control. It determines how quickly the circuit is changing over time (by looking at the derivative of the error signal) and multiplies it by K_d to produce the derivative response.

$$D = K_d \frac{d}{dt}e(t)$$

Unlike proportional and integral control, derivative control will slow the response of the circuit. In doing so, it is able to partially compensate for the overshoot as well as damp out any oscillations caused by integral and

proportional control. High gain values cause the circuit to respond very slowly and can leave one susceptible to noise and high frequency oscillation (as the circuit becomes too slow to respond quickly). Too low and the circuit is prone to overshooting the set point value. However, in some cases overshooting the set point value by any significant amount must be avoided and thus a higher derivative gain (along with lower proportional gain) can be used. The chart below explains the effects of increasing the gain of any one of the parameters independently.

Parameter Increased	Rise Time	Overshoot	Settling Time	Steady-State Error	Stability
K_p	Decrease	Increase	Small Change	Decrease	Degrade
K_i	Decrease	Increase	Increase	Decrease Significantly	Degrade
K_d	Minor Decrease	Minor Decrease	Minor Decrease	No Effect	Improve (for small K_d)

5.3 Bumpless Integration

When the integration term K_i is applied outside the integral ($I = K_i \int_0^t e(\tau)d\tau$) and K_i is zero, the value of the integral can accumulate to very large values if *e(t)* maintains the same sign for an extended period of time. When K_i is then set to nonzero, this can result in a sudden "bump" in the output value as the accumulated value of the integral is suddenly able to affect the output.

The DSC1 controller instead applies the integration term inside the integral ($I = \int_0^t K_i(\tau)e(\tau)d\tau$). In this way when K_i is zero the value of the integral stays at zero and there is no bump (thus "bumpless") when K_i is set to a nonzero value.

Notice the DSC1 uses a value of zero for K_i when the I term is disabled, regardless of the value entered by the user.

5.4 Discrete-Time Servo Controllers

5.4.1 Data Format

The "error" signal, $e(t)$, in the PID loop at time t is determined as

$$e(t) = S - v(t)$$

Where S is the setpoint voltage and $v(t)$ is the voltage sample at time t.

5.4.2 Control Law in the Time Domain

Three gain terms are computed and summed together.

$$u(t) = G_p(t) + G_i(t) + G_d(t)$$

$$G_p(t) = K_p e(t)$$

$$G_i(t) = \int_0^t K_i(\tau)e(\tau)d\tau$$

$$G_d(t) = K_d \frac{d}{dt} e(t)$$

Where $G_p(t)$, $G_i(t)$, and $G_d(t)$ are the proportional, integral, and derivative gains, respectively. K_p, K_i, and K_d are the proportional, integral, and derivative gain coefficients.

5.4.3 Approximating the Integral and the Derivative

The DSC1 controller approximates an integrator with an accumulator.

$$\int_0^t K_i(\tau)e(\tau)d\tau \approx h\sum_{n=1}^{N} K_i\,[n]e[n]$$

Where *n* is the index of the point acquired at time *t* and h is the time between ADC samples.

We make a similar approximation to the derivative using the difference between $e[n]$ and $e[n-1]$.

$$\frac{d}{dt}e(t) \approx \frac{1}{h}(e[n] - e[n-1])$$

Therefore, we adjust our prior derivation of the control law to:

$$u[n] = K_p e[n] + h\sum_{n=1}^{N} K_i\,[n]e[n] + \frac{K_d}{h}(e[n] - e[n-1])$$

5.5 Servo Mode: PID Tuning

In general, the gains of P, I, and D will need to be adjusted by the user in order to optimize the performance of the system. While there is not a static set of rules for what the values should be for any specific system, following the general procedures should help in tuning a circuit to match one's system and environment. In general, a properly tuned PID circuit will typically overshoot the setpoint value slightly and then quickly damp out to reach the setpoint value and hold steady at that point.

The PID loop can lock to either a positive or negative slope by changing the sign of the P, I, and D gains. In the DSC1, the signs are locked together so changing one will change them all.

Manual tuning of the gain settings is the simplest method for setting the PID controls. However, this procedure is done actively (the PID controller attached to the system and the PID loop enabled) and requires some amount of experience to achieve good results.

To tune your PID controller manually, first set the integral and derivative gains to zero. Increase the proportional gain until you observe oscillation in the output. Your proportional gain should then be set to roughly half this value. After the proportional gain is set, increase the integral gain until any offset is corrected for on a time scale appropriate for your system. If you increase this gain too much, you will observe significant overshoot of the SP value and instability in the circuit. Once the integral gain is set, the derivative gain can then be increased. Derivative gain will reduce overshoot and damp the system quickly to the set point value. If you increase the derivative gain too much, you will see large overshoot (due to the circuit being too slow to respond). By playing with the gain settings, you can optimize the performance of your PID circuit, resulting in a system that quickly responds to changes and effectively damps out oscillation about the set point value.

Control Type	K_p	K_i	K_d
P	0.50 K_u	-	-
PI	0.45 K_u	1.2 K_p/P_u	-
PID	0.60 K_u	2 K_p/P_u	$K_pP_u/8$

While manual tuning can be very effective at setting a PID circuit for your specific system, it does require some amount of experience and understanding of PID circuits and response. The Ziegler-Nichols method for PID tuning offers a more structured guide to setting PID values. Again, you'll want to set the integral and derivative gain to zero. Increase the proportional gain until the circuit starts to oscillate. We will call this gain level K_u. The oscillation will have a period of P_u. Gains are for various control circuits are then given in the chart above.

Note that when using the Ziegler-Nichols tuning method with the DSC1, the integral term determined from the table should be multiplied by $2 \cdot 10^{-6}$ to normalize to the sample rate. Similarly, the derivative coefficient should be divided by $2 \cdot 10^{-6}$ to normalize to the sample rate.

5.6 Ramp Mode

Users may often need to determine the large-signal operating point or useful setpoint for a system. To determine either the large-signal operating point (hereafter referred to as DC offset) or optimal servo setpoint, a common technique is to simply stimulate the system repeatedly with a linearly increasing voltage signal. The pattern is commonly referred to as a sawtooth wave, for its resemblance to the teeth of a saw. After each ramp the control variable is returned to its minimum value to allow the system to settle before the start of the next ramp cycle.

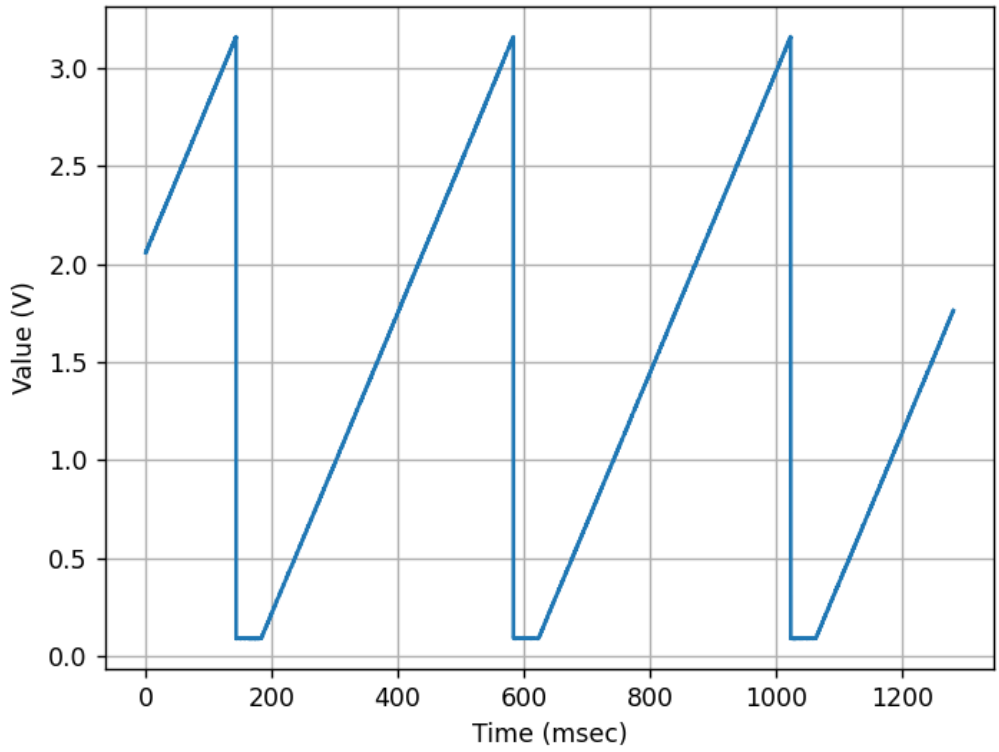


Figure 12 Oscillogram of a Several Cycles of the Ramp Mode Waveform with Amplitude Set to 1.5 V and Offset Set to 1.65 V

5.7 Peak Lock Mode

The peak lock mode implements a dither locking algorithm also known as an extremum seeking controller. In this mode of operation, the control value is superimposed on a sine wave output. The measured input voltage is first digitally high-pass filtered (HPF) to remove any DC offset. Then the AC coupled signal is demodulated by multiplying each measured voltage by the outgoing sine wave modulation value. This multiplication operation creates a demodulated signal with two principal components: a sine wave at the sum of the two frequencies and a signal at the difference of the two frequencies. A second digital filter, this time a low pass filter (LPF), attenuates

the sum-of-two frequencies signal, and transmits the low frequency difference-of-two frequencies signal. Signal content at the same frequency as the modulation appears as a DC signal post demodulation. The final step in the peak lock algorithm is to integrate the LPF signal. The integrator output, combined with the outgoing modulation, drives the output voltage. The accumulation of low frequency demodulated signal energy in the integrator pushes the offset control variable of the output higher and higher until the sign of the LPF output reverses and the integrator output begins to diminish. As the control value approaches the peak of the system response, the result of the modulation on the input signal to the servo controller becomes smaller and smaller, since the slope of a sinusoidal wave form is zero at its peak. This in turn means that there is a lower output value from the low-pass-filtered, demodulated signal, and therefore less to accumulate in the integrator.

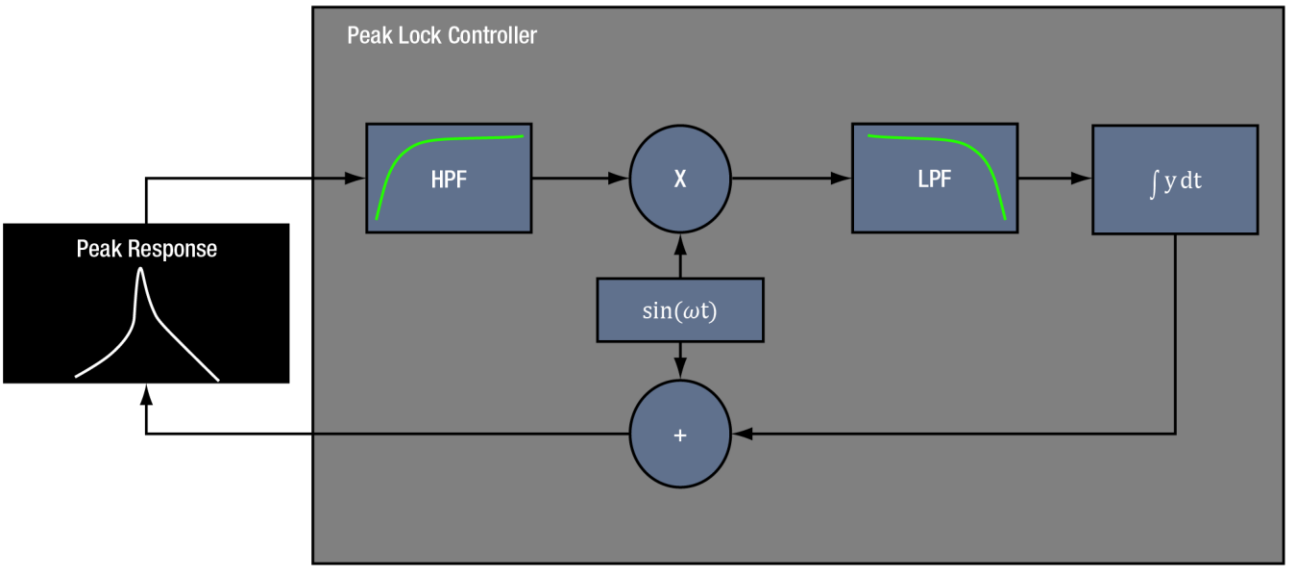


Figure 13 Block diagram of a peak locking controller. The input signal from the plant is digitized, then high-pass filtered. The HPF output signal is demodulated with a digital local oscillator. The output of the demodulator is low-pass filtered and then integrated. The integrator output is added to the modulation signal and output to the plant.

Peak locking is a good control algorithm to choose when the system the user wishes to control does not have a monotonic response around the optimal control point. Examples of these sorts of systems are optical media with a resonant wavelength, such as a vapor cell, or an RF band-reject filter (notch filter). The central characteristic of the peak locking control scheme is the algorithm's tendency to steer the system towards the zero-crossing of the error signal which coincides with a peak in the measured signal, as if the error signal were the derivative of measured signal. Note that the peak may be positive or negative.

To get started with the peak locking mode of operation for the DSC1, you may follow this procedure.

1. Make sure that there is a peak (or valley) of the signal you are locking to is within the control variable range of the actuator, and that the peak position is relatively stable with time. It is helpful to use the RAMP mode to visualize the signal over the control variable range of interest.
2. Note the control variable position of the peak (or valley).

3. Estimate how broad the peak (or valley) is in control variable at half of the height of the peak. This width, in volts, is commonly referred to as the Full-Width Half-Max or FWHM. It should be at least 0.1V wide for good results.
4. Set the modulation amplitude (A) to 1% to 10% of the FWHM voltage.
5. Set the offset voltage as close as possible to the position of the peak (or valley) that you desire to lock to.
6. Set the modulation frequency to the desired frequency. On the touch screen this is affected through the M, modulation frequency parameter. The modulation frequency is 100 Hz times M. The best modulation frequency selection depends on the application. Thorlabs recommends values around 1 kHz for mechanical actuators. Higher frequencies may be used to electro-optic actuators.
7. Set the peak lock integral coefficient (K) to 0.1 times A. K can be positive or negative. Generally, positive K locks to a peak of the input signal, while negative K locks to a valley of the input signal. However, if the actuator or system being locked has more than 90 degree phase delay at the dither frequency, the sign of K will invert and positive K will lock to a valley, and negative K will lock to a peak.
8. Press Run and verify that the control variable output changes from the original offset (O) value and does not run away to an extreme. Alternatively, monitor the process variable using an oscilloscope to verify that the DSC1 is locking to the peak or valley desired.

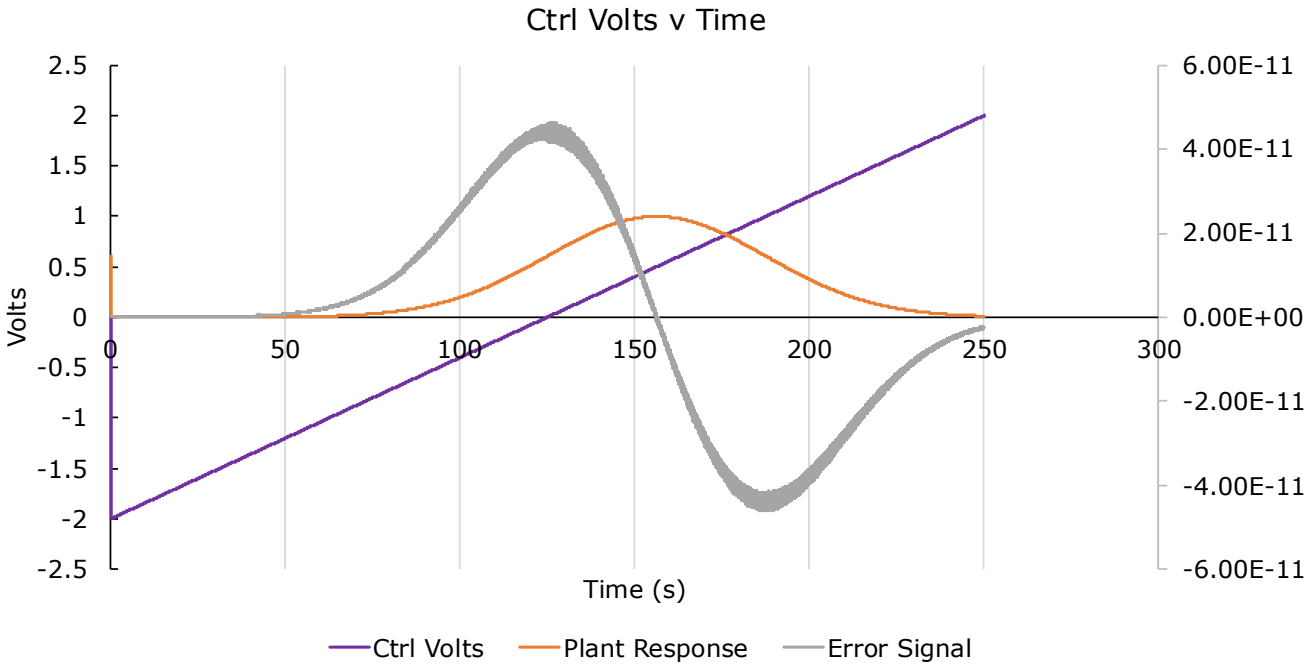


Figure 14 Example data from ramping the output offset voltage over a plant, with response shown in orange. The error signal is the calculated value after the low pass filter in Figure 13. Note the error signal zero crossing aligns with the peak of the plant response signal, which will be the lock point for PEAK mode.

Chapter 6 Maintenance and Cleaning

The DSC1 requires no regular maintenance. Should the touchscreen on the device become dirty, Thorlabs recommends gently cleaning the touchscreen with a soft, lint-free cloth, saturated with diluted isopropyl alcohol.

Chapter 7 Troubleshooting and Repair

The table below describes typical issues with the DSC1 and Thorlabs recommended remedies.

Issue	Explanation	Remedy
The device does not turn on when plugged into USB Type-C power.	The device requires as much as 750 mA of current from a 5 V supply, 3.75 W. This may exceed the power capabilities of some USB-A connectors on laptops and PCs.	Use Thorlabs DS5 or CPS1 power supplies. Alternatively, use a USB Type-C power supply such as is typically used to charge a phone or laptop that is rated to output at least 750 mA at 5 V.
The device does not turn on when the data port is plugged into a PC.	The DSC1 only draws power from the USB Type-C power connector. The USB Type Mini-B connector is data only.	Connect the USB Type-C port to a power supply rated to output at least 750 mA at 5 V, such as Thorlabs DS5 or CPS1.
Peak of the signal shown in Ramp mode does not correspond exactly with the lock point in Peak mode.	When modulating a system (like a laser diode current), there are often multiple time scales for the output signal to affect the device. Changing the current will impart both a fast (current density) and slow (temperature) shift in the laser wavelength. This can lead the location of the peak in the Ramp scan to be dependent on the average scan voltage (offset) as well as the instantaneous voltage.	Adjust the system or laser to center the location of the peaked feature at the average value of the scan. This should minimize the difference between the Ramp peak voltage and the steady-state locked voltage.

Chapter 8 Disposal



Thorlabs verifies our compliance with the WEEE (Waste Electrical and Electronic Equipment) directive of the European Community and the corresponding national laws. Accordingly, all end users in the EC may return “end of life” Annex I category electrical and electronic equipment sold after August 13, 2005 to Thorlabs, without incurring disposal charges. Eligible units are marked with the crossed out “wheelie bin” logo (see right), were sold to and are currently owned by a company or institute within the EC and are not dissembled or contaminated. Contact Thorlabs for more information. Waste treatment is your own responsibility. “End of life” units must be returned to Thorlabs or handed to a company specializing in waste recovery. Do not dispose of the unit in a litter bin or at a public waste disposal site. It is the user’s responsibility to delete all private data stored on the device prior to disposal.

Chapter 9 Thorlabs Worldwide Contacts

For technical support or sales inquiries, please visit us at www.thorlabs.com/contact for our most up-to-date contact information.



Corporate Headquarters
Thorlabs, Inc.
43 Sparta Ave
Newton, New Jersey 07860
United States

sales@thorlabs.com
techsupport@thorlabs.com

Product Manufacturer
Thorlabs, Inc.
43 Sparta Ave
Newton, New Jersey 07860
United States

sales@thorlabs.com
techsupport@thorlabs.com

EU Importer
Thorlabs GmbH
Münchner Weg 1
D-85232 Bergkirchen
Germany

sales.de@thorlabs.com
europe@thorlabs.com

UK Importer
Thorlabs Ltd.
204 Lancaster Way Business Park
Ely CB6 3NX
United Kingdom

sales.uk@thorlabs.com
techsupport.uk@thorlabs.com



THORLABS

www.thorlabs.com