PDAPC7
Extended InGaAs Fixed Gain Photodetector, OEM Package

User Guide



Table of Contents

Chapter 1 Introduction 1

1.1 Intended Use 1

1.2 Explanation of Safety Warnings (Required if Used in Manual) 1

1.3 Description 1

1.4 Technical Data 2

1.4.1 Specifications 2

1.4.2 Graphs 3

1.4.3 Mechanical Drawing 4

1.5 Pin Diagrams 4

Chapter 2 Installation 5

2.1 Installation Instructions 5

Chapter 3 Operation 7

3.1 Theory of Operation 7

3.2 Responsivity 7

3.3 Dark Current 7

3.4 Bandwidth and Response 8

3.5 Terminating resistance 8

Chapter 4 Troubleshooting and Repair 9

Chapter 5 Disposal 9

Chapter 6 Thorlabs Worldwide Contacts (Required) 10

Chapter 1 Introduction

1.1 Intended Use

This product is intended to be used to monitor near infrared (NIR) light levels and create an appropriate electrical signal in a laboratory environment. The OEM style PCB is intended to be used in tight places, with the Thorlabs 30mm cage systems, or embedded in user systems.

The photodiode is centered on the 30mm spacing. These mounting holes can be used for mounting into custom housings and devices. A 10-pin header is provided to supply power and retrieve the detector voltage.

1.2 Explanation of Safety Warnings

NOTICE Indicates information considered important, but not hazard-related, such as possible damage to the product.


Shock Warning


ESD Component Caution


The wheelie bin symbol on the product, the accessories or packaging indicates that this device must not be treated as unsorted municipal waste but must be collected separately.

1.3 Description

The PDAPC7 is a fixed-gain, Extended Indium Gallium Arsenide (InGaAs) photodetector in an OEM package with an operating wavelength range of 900 to 2600 nm. A buffered output drives 50 Ω load impedances up to 5 V. The output signal is readable from the MMCX connector. MMCX to BNC cables (Item # CA3339) or MMCX to SMA cables (Item # CA3439) are available for purchase separately. Users can also access the digital pins of the output through an on-board header.

NOTICE This product is ESD sensitive. Take all appropriate precautions to discharge personnel and equipment before making any electrical connections to the unit.

1.4 Technical Data

1.4.1 Specifications


Electrical Specifications[a]		
Detector	-	Extended InGaAs
Active Area	-	Ø1.0 mm (0.8 mm^2)
Wavelength Range	λ	900 - 2600 nm
Peak Wavelength (Typical)	λ_p	2300 nm
Peak Response (Typical)	$\Re(\lambda_p)$	1.35 A/W
Small Signal Bandwidth	-	25 MHz
Rise Time	-	15 ns
Noise Equivalent Power[b]	-	1.01 x 10^{-11} W/√Hz
Noise (RMS)	-	2.5 mV
Max Output Current	I_{OUT}	100 mA
Dark Offset	-	±75 mV (Typ.) ±375 mV (Max)
Load Impedance	-	50 Ω to Hi-Z
Transimpedance Gain	-	5 x 10^3 ± 2% V/A (50 Ω) 1 x 10^4 ± 2% V/A (Hi-Z)
Output Voltage	V_{OUT}	0 to 5 V (50 Ω) 0 to 10 V (Hi-Z)

a. The PDAPC7 has a 50 Ω series terminator resistor (i.e., in series with the amplifier output). This forms a voltage divider with any load impedance (e.g., 50 Ω load divides signal in half).

b. The Noise Equivalent Power is specified at the peak responsivity wavelength.

General	
PCB Size	1.47" x 1.47" (37.3 mm x 37.3 mm)
Operating Temperature	10 to 40 °C
Storage Temperature	-20 to 70 °C

1.4.2 Graphs

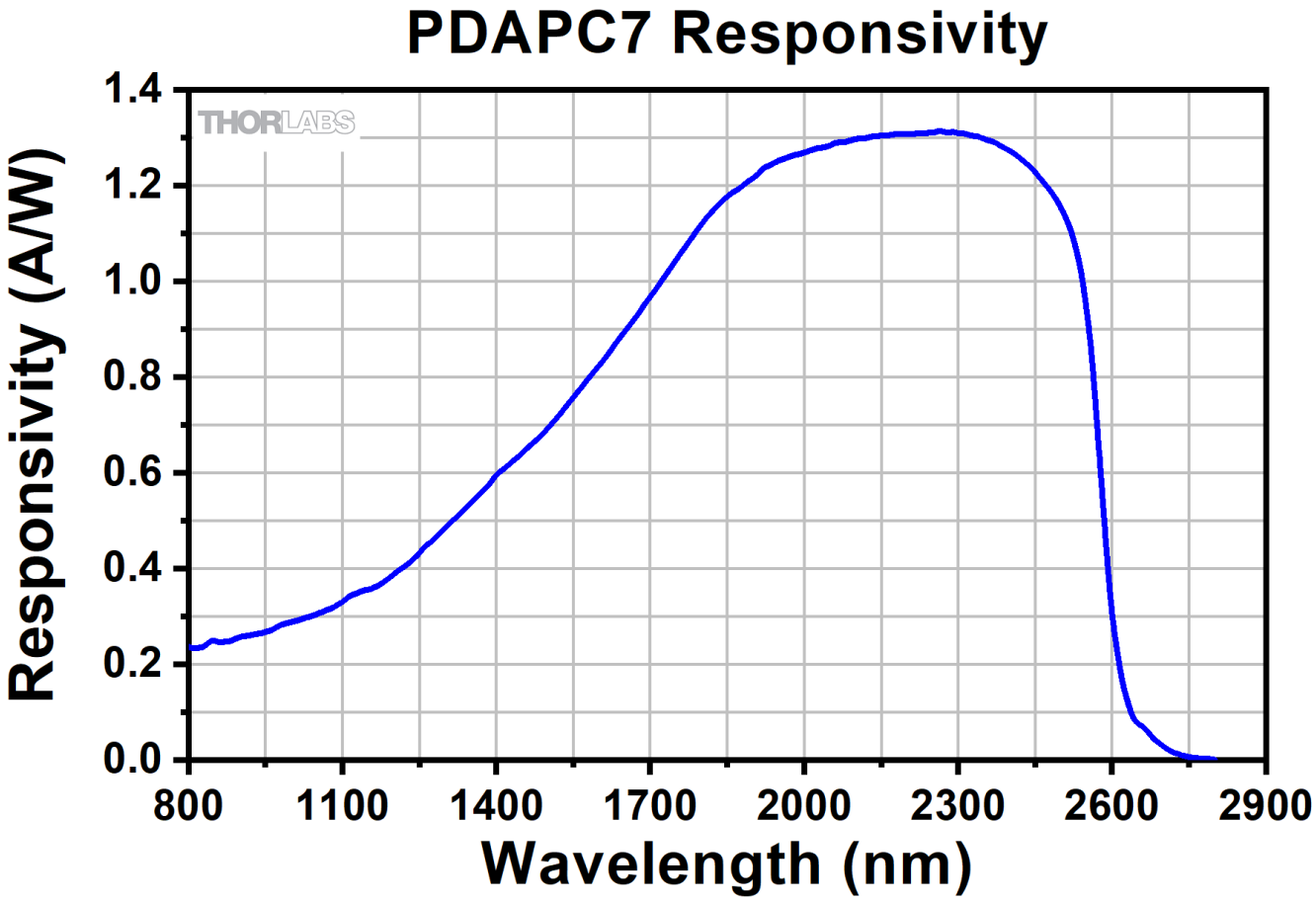
PDAPC7 Responsivity
THORLABS
Responsivity (A/W)
1.4
1.2
1.0
0.8
0.6
0.4
0.2
0.0
800
1100
1400
1700
2000
2300
2600
2900
Wavelength (nm)

1.4.3 Mechanical Drawing

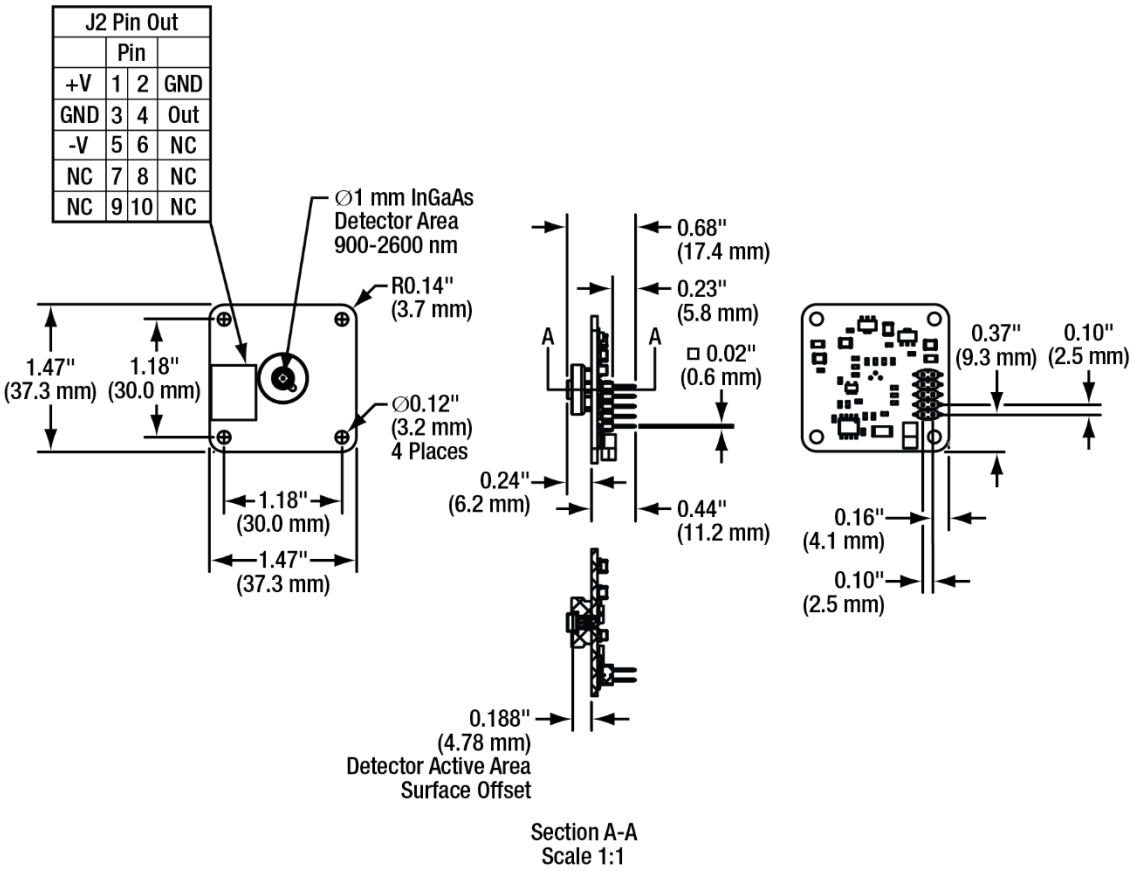


Figure 1 Mechanical Drawing

1.5 Pin Diagrams

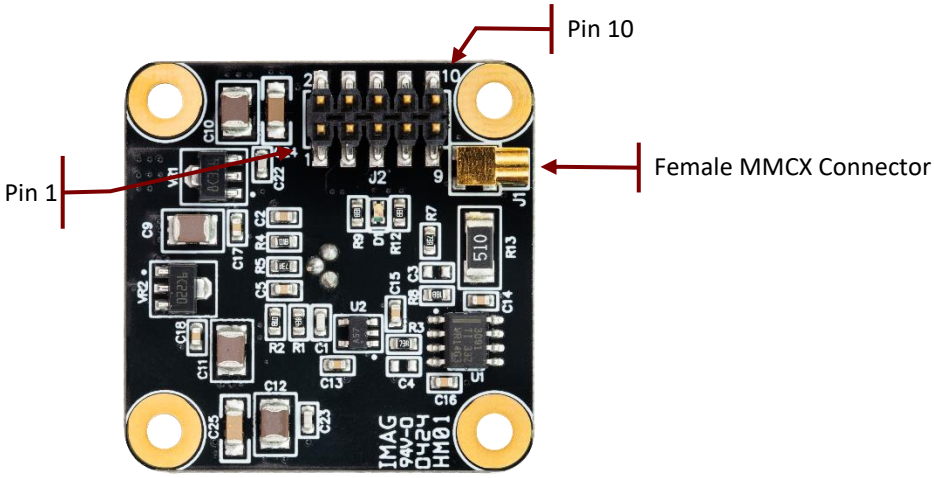


Figure 2 J2 Header

1. The pin out table for the on-board header (J2) is shown in the table below. Connect a power supply to pins 1, 3, and 5 of the J2 header. The required power supply voltages are:

J2 Pin Out[a]			
Pin	**Designation**	**Pin**	**Designation**
1[a]	+12 VDC	**2**	GND
3	GND	**4**	OUTPUT
5[a]	-12 VDC	**6**[b]	NC
7[b]	NC	**8**[b]	NC
9[b]	NC	**10**[b]	NC

a. ±V refers to the power supply positive/negative voltages.

b. No Connection; Reserved for Switchable Gain Devices

Chapter 2 Installation

2.1 Installation Instructions

1. Connect a power supply to pins 1, 3, and 5 of the J2 header. The required power supply voltages are:

J2 Pin	Voltage
1	+12 VDC
3	GND
5	-12 VDC

Thorlabs offers the LD1255-SUPPLY as one compatible power supply option. Connection to the board is shown below:

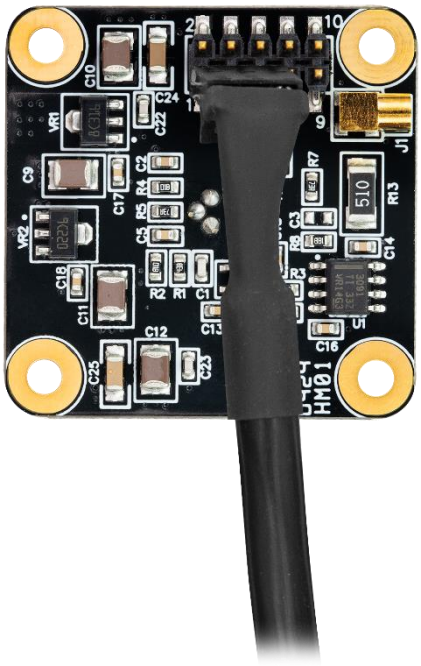


Figure 3 Power Supply Connection

2. Attach a 50 Ω coax cable to the MMCX output of the PDAPC7 detector. Thorlabs offers two suitable cable options, the CA3339 MMCX to BNC cable or the CA3439 MMCX to SMA cable. When running cable lengths longer than 12”, Thorlabs recommends terminating the opposite end of the coax with a 50 Ω resistor for optimal performance. Thorlabs Item # T4119 is suitable for this purpose. Connect the remaining end to a measurement devise such as an oscilloscope or high-speed DAQ card.


NOTICE Many high speed oscilloscopes have input impedances of 50 Ω. In this case, do not install a 50 Ω terminator. The combined loads will equal 25 Ω which could allow ~135 mA of output current. This will damage the output driver of the PDA.

3. Power the PDAPC7 on. When the unit is on, the power indicator LED will be illuminated yellow. If the LED is green or red, this indicates one power rail, either positive or negative, respectively, is connected. In that case, the user would need to check the power inputs to make sure both power rails are on and connected correctly.
4. Apply a light source to the detector.
5. During alignment, the user should take appropriate precautions. For example, ensuring the incident power is within safe operating power. Additionally, the user should have proper eye and skin protection as recommended by the light source manufacturer.

Chapter 3 Operation

3.1 Theory of Operation

Thorlabs PDA series are ideal for measuring both pulsed and CW light sources. The PDAPC7 includes a reverse-biased PIN photodiode, mated to a fixed gain transimpedance amplifier.

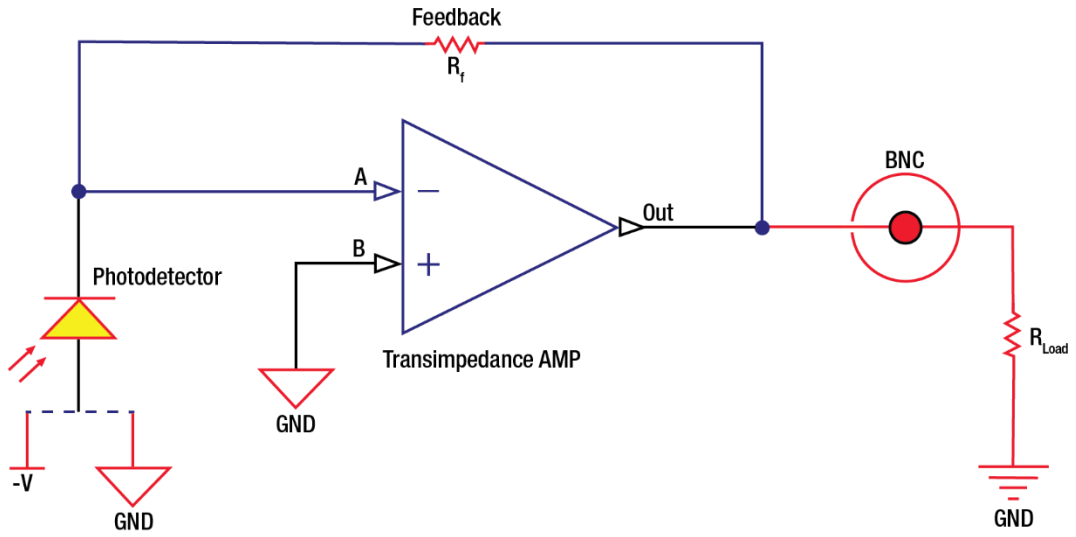

3.2 Responsivity

The responsivity of a photodiode can be defined as a ratio of generated photocurrent (IPD) to the incident light power (P) at a given wavelength:

$$R(\lambda) = \frac{I_{PD}}{P}$$

3.3 Dark Current

Dark current is leakage current which flows when a bias voltage is applied to a photodiode. The PDA with Transimpedance Amplifier does control the dark current flowing out. Looking at the figure above, it can be noted that Point B is held at ground and the amplifier will try to hold point A to “Virtual Ground”. This minimizes the effects of dark current present in the system.

The dark current present is also affected by the photodiode material and the size of the active area. Silicon devices generally produce low dark current compared to germanium devices which have high dark currents. The table below lists several photodiode materials and their relative dark currents, speeds, sensitivity, and costs. Please note that sensitivity values in the table are typical values; Thorlabs offers photodetectors with sensitivity ranges that vary from those shown below.

Material	Dark Current	Speed	Sensitivity (nm)[1]	Cost
Silicon (Si)	Low	High	400 - 1000	Low
Germanium (Ge)	High	Low	900 - 1600	Low
Gallium Phosphide (GaP)	Low	High	150 - 550	Med
Indium Gallium Arsenide (InGaAs)	Low	High	800 - 1800	Med
Extended Range: Indium Gallium Arsenide (InGaAs)	High	High	1200 - 2600	High
Lead Sulfide (PbS)	N/A[2]	Low	1000 - 2900	Med
Lead Selenide (PbSe)	N/A[3]	Low	1500 - 4800	Med
Indium Arsenide Antimonide (InAsSb)	N/A[4]	High	1000 - 5300[5]	Med

3.4 Bandwidth and Response

For the PDA detectors, the gain of the detector is dependent on the feedback element (R_f). The bandwidth of the detector can be calculated using the following:

$$f(-3dB) = \sqrt{\frac{GBP}{4\pi R_f \times C_D}}$$

where GBP is the amplifier gain bandwidth product and C_D is the sum of the photodiode junction capacitance and the amplifier capacitance.

3.5 Terminating Resistance

The maximum output of the PDAPC7 is 10 V for high impedance loads (i.e. R_{Load} > 5 kΩ) and 5 V for 50 Ω loads. Adjust the gain so that the measured signal level out of the PDAPC8 is below 10 V (5 V with a 50 Ω load) to avoid saturation.

For low terminating resistors, <5 kΩ or 1% error, an additional factor needs to be considered. The output of the PDA includes a 50 Ω series resistor (R_s). The output load creates a voltage divider with the 50 Ω series resistor as follows:e

$$Scale\ Factor\ = \frac{R_{Load}}{R_{Load} + R_s}$$

$$V_{OUT}\ =\ \mathcal{R}(\lambda) * Transimpedance\ Gain * Scale\ Factor * Input\ Power\ (W)$$

Note that we already include the scale factor in our specification for the gain at 50 Ω. Refer to the table in Chapter 6 for additional performance specifications.

[1] *Approximate values; actual wavelength values will vary.*
[2] *These are photoconductive detectors, normally AC coupled.*
[3] *These are photoconductive detectors, normally AC coupled.*
[4] *Photovoltaic detectors, generally unbiased.*
[5] *At 25 °C.*

Chapter 4 Troubleshooting and Repair

Problem	Suggested Solutions
There is no signal response.	Verify that the power is switched on and all connections are secure.
	Verify the proper terminating resistor is installed if using a Voltage measurement device.
	Verify that the optical signal wavelength is within the specified wavelength range.
	Verify that the optical signal is illuminating the detector active area.
Output voltage will not increase. **Detector output is skewed.**	Check to make sure the detector is not saturated. Refer to the Output Voltage specification in the table on page 2.

Chapter 5 Disposal



Thorlabs verifies our compliance with the WEEE (Waste Electrical and Electronic Equipment) directive of the European Community and the corresponding national laws. Accordingly, all end users in the EC may return "end of life" Annex I category electrical and electronic equipment sold after August 13, 2005 to Thorlabs, without incurring disposal charges. Eligible units are marked with the crossed out "wheelie bin" logo (see right), were sold to and are currently owned by a company or institute within the EC and are not dissembled or contaminated. Contact Thorlabs for more information. Waste treatment is your own responsibility. "End of life" units must be returned to Thorlabs or handed to a company specializing in waste recovery. Do not dispose of the unit in a litter bin or at a public waste disposal site. It is the user's responsibility to delete all private data stored on the device prior to disposal.

Chapter 6 Thorlabs Worldwide Contacts

For technical support or sales inquiries, please visit us at www.thorlabs.com/contact for our most up-to-date contact information.



Corporate Headquarters
Thorlabs, Inc.
43 Sparta Ave
Newton, New Jersey 07860
United States

sales@thorlabs.com
techsupport@thorlabs.com

Product Manufacturer
Thorlabs, Inc.
43 Sparta Ave
Newton, New Jersey 07860
United States

sales@thorlabs.com
techsupport@thorlabs.com

EU Importer
Thorlabs GmbH
Münchner Weg 1
D-85232 Bergkirchen
Germany

sales.de@thorlabs.com
europe@thorlabs.com

UK Importer
Thorlabs Ltd.
204 Lancaster Way Business Park
Ely CB6 3NX
United Kingdom

sales.uk@thorlabs.com
techsupport.uk@thorlabs.com



THORLABS

www.thorlabs.com