Slim Dual Range Power Head with Silicon Detector



S130VC

Description

The *S130VC* with its slim design of only 5 mm at the detector side is perfect to fit into tight optical setups. The ruggedized aluminum head with its large active area silicon photodiode can be held by hand or can be mounted with its 8-32 and M4 threaded mounting holes to posts and post holders.

A slideable ND filter enlarges the power range by approximately a factor of 100; the position of the filter is automatically recognized by the connected console after a few seconds for calculating the power level with the right responsivity value.

For sensor readout, the S120VC is compatible with any Thorlabs consoles and interfaces equipped with a DB9 connection, and rated for photodiode power measurement. A non-volatile memory in the sensor connector contains sensor information data and the NIST and PTB traceable calibration data.

Available Accessories

SM1A29	SM1 Thread Adapter with Alignement Target
S120-FC	FC Fiber Adapter
S120-SMA	SMA Fiber Adapter
S120-SC	SC Fiber Adapter
S120-LC	LC Fiber Adpter
S120-ST	ST Fiber Adapter

The *S130VC* is also compatible to the Thorlabs imperial and metric post and post-holder series and with the optional thread adapter to Thorlabs SM1 mechanics.

Cleaning and Maintenance

There are no serviceable parts in the *S130VC* head. The housing may be cleaned by wiping with a soft damp cloth. When cleaning the aperture filter, treat it as any other fine optic. Gently blow off any debris using compressed air and wipe gently with an optic tissue wetted with propanol. If you suspect a problem with your *S130VC* please call Thorlabs and an engineer will be happy to assist you.

As long as the sensor has not been exposed to excessive optical power (please pay attention to the maximum ratings in the technical specifications), the calibration should be very stable over long periods of time (well over a year). To keep the accuracy and performance of the *S130VC*, Thorlabs recommends a yearly recalibration, starting one year after purchase.


Specifications

Specification	Value
Detector Type	Silicon Photodiode (UV extended)
Wavelength Range	200 - 1100 nm
Optical Power Working Range	500 pW - 0.5 mW (Up to 50 mW) (200 - 450 nm) 500 pW - 0.2 mW (Up to 20 mW) (450 – 1100 nm)
Max Average Power Density	20 W/cm^2
Max Pulse Energy	20 µJ
Linearity	± 0.5%
Resolution[a]	100 pW
Measurement Uncertainty[b]	±3% (440 - 980 nm) ±5% (280 - 439 nm) ±7% (200 - 279 nm, 981 - 1100 nm)
Typical Application	Low Power Lasers
Laser Types	Diode, Diode Arrays, He-Ne, Dye, Ion Lasers (Ar+, Kr+)
Coating /Diffuser[c]	Reflective ND (OD1.5)
Cooling	Convection
Response Time	< 1 µs
Sensor Dimensions	150 mm x 19 mm x 10 mm 5 mm Thickness on Sensor Side
Active Detector Area	9.7 mm x 9.7 mm
Input Aperture	Ø9.5 mm
Cable Length	1.5 m
Connector	Sub-D 9p Male
Weight	0.125 kg
Post	#8-32 & M4 Threads
Aperture Thread (Optional)	SM1, Outer Thread with SM1A29 Adapter
Fiber Adapters (Optional)	FC, SC, LC, SMA, ST (SM1A29 Adapter Required)


a. Measured with PM100D console in bandwidth low setting without filter.
b. Beam Diameter > 1mm
c. This specification is valid for S130VC devices from serial number 1203xxx. For older versions, please contact technical support.

Please note that the S130VC power meter head is not compatible with the older Thorlabs power meter consoles (PM100, PM30, PM300, PM300E, S100).

Typical Response Graph

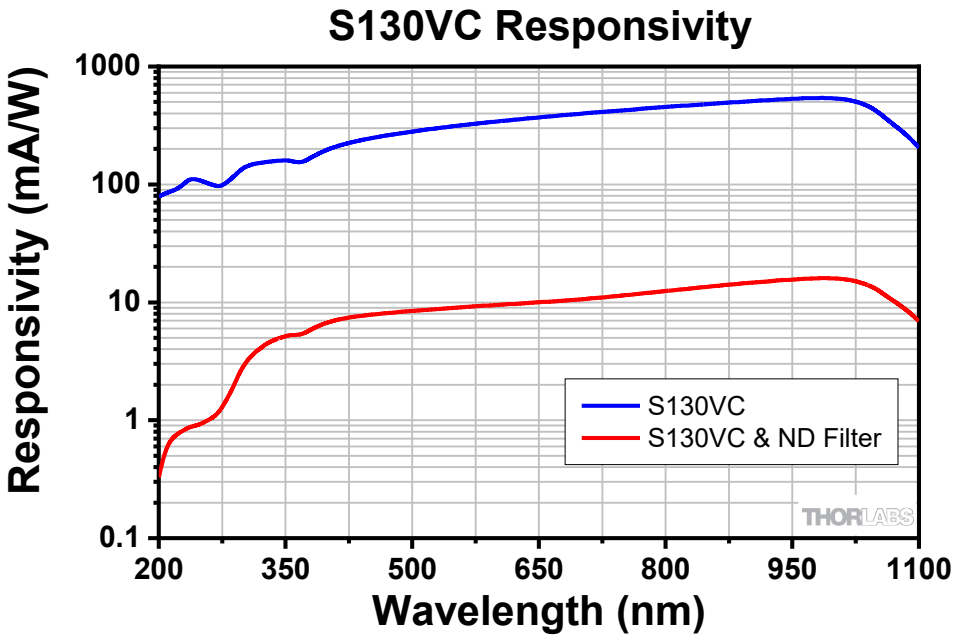



Drawings

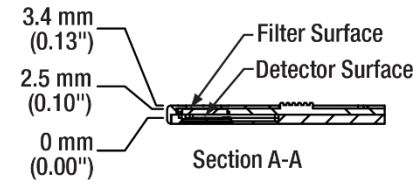

Connector Pin Diagram

2	EEPROM Data
3	PD Anode
4	PD Kathode
6	EEPROM Ground
7	Slider Detection
1,4,5,9	N.C.

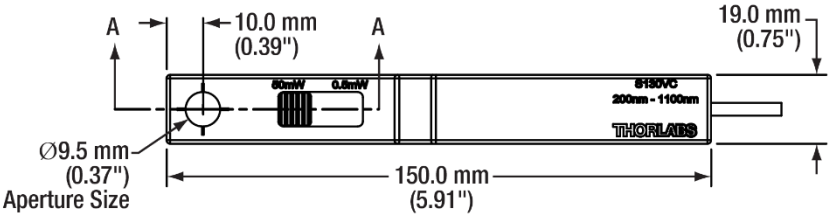

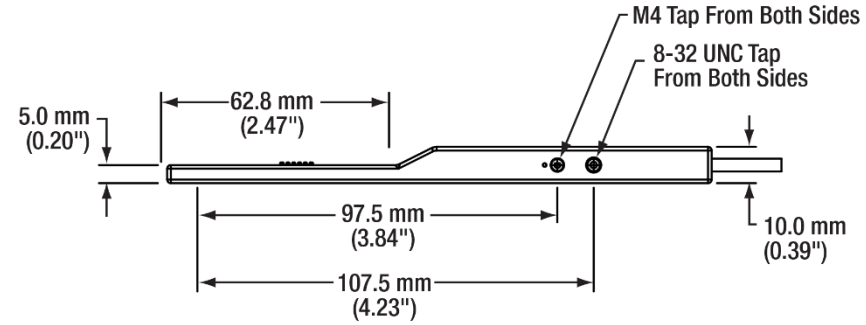

Precautions and Warranty Information

These products are covered under a two year warranty. In order to ensure the proper functioning of a sensor care must be given to maintain the highest standards of compliance to the maximum optical specifications when handling such devices. The photodiodes are particularly sensitive to any value that exceeds the absolute maximum ratings of the product. Any applied optical power in excess of the maximum specification will cause damage and possible complete failure to the product.

Thorlabs, Inc. Life Support and Military Use Application Policy is stated below:

THORLABS' PRODUCTS ARE NOT AUTHORIZED FOR USE AS CRITICAL COMPONENTS IN LIFE SUPPORT DEVICES OR SYSTEMS OR IN ANY MILITARY APPLICATION WITHOUT THE EXPRESS WRITTEN APPROVAL OF THE PRESIDENT OF THORLABS, INC. As used herein:

1. Life support devices or systems are devices or systems which, (a) are intended for surgical implant into the body, or (b) support or sustain life, and whose failure to perform, when properly used in accordance with instructions for use provided in the labeling, can be reasonably expected to result in a significant injury to the user.

2. A critical component is any component in a life support device or system whose failure to perform can be reasonably expected to cause the failure of the life support device or system or to affect its safety or effectiveness.

3. The Thorlabs products described in this document are not intended nor warranted for usage in Military Applications.

